### Ordering Details



Iltimate Endoscopic Linear Cutter Staplers with Universal Handles								
Product	Code	Shaft Length	Pack QTY					
<b>Ultimate</b> Endoscopic Linear Cutter Stapler with Universal Handle	PS9ECHU06	60mm	3 / Box					
Ultimate Endoscopic Linear Cutter Stapler with Universal Handle	PS9ECHU16	160mm	3 / Box					
<b>Ultimate</b> Endoscopic Linear Cutter Stapler with Universal Handle	PS9ECHU26	260mm	3 / Box					

Ultimate Endoscopic Linear Cutter Stapler Articulating Reloading Units											
Product	Code	Colour	Tissue Type	Cut Length	Rows	Staples	Open Staple Height	Closed Staple Height	Pack QTY		
<b>Ultimate</b> 30mm Endoscopic Articulating Reloading Unit	PS9ECAR30W	White	Thin	25mm	6	48	2.5mm	1.0mm	6 / Box		
Ultimate 30mm Endoscopic Articulating Reloading Unit	PS9ECAR30B	Blue	Standard	25mm	6	48	3.5mm	1.5mm	6 / Box		
Ultimate 45mm Endoscopic Articulating Reloading Unit	PS9ECAR45W	White	Thin	40mm	6	66	2.5mm	1.0mm	6 / Box		
<b>Ultimate</b> 45mm Endoscopic Articulating Reloading Unit	PS9ECAR45B	Blue	Standard	40mm	6	66	3.5mm	1.5mm	6 / Box		
<b>Ultimate</b> 45mm Endoscopic Articulating Reloading Unit	PS9ECAR45D	Gold	Regular - Thick	40mm	6	66	4.2mm	1.8mm	6 / Box		
<b>Ultimate</b> 45mm Endoscopic Articulating Reloading Unit	PS9ECAR45G	Green	Thick	40mm	6	66	4.8mm	2.0mm	6 / Box		
<b>Ultimate</b> 45mm Endoscopic Articulating Reloading Unit	PS9ECAR45N	Black	Very Thick	40mm	6	66	5.0mm	2.35mm	6 / Box		
Ultimate 60mm Endoscopic Articulating Reloading Unit	PS9ECAR60W	White	Thin	55mm	6	90	2.5mm	1.0mm	6 / Box		
Ultimate 60mm Endoscopic Articulating Reloading Unit	PS9ECAR60B	Blue	Standard	55mm	6	90	3.5mm	1.5mm	6 / Box		
Ultimate 60mm Endoscopic Articulating Reloading Unit	PS9ECAR60D	Gold	Regular - Thick	55mm	6	90	4.2mm	1.8mm	6 / Box		
Ultimate 60mm Endoscopic Articulating Reloading Unit	PS9ECAR60G	Green	Thick	55mm	6	90	4.8mm	2.0mm	6 / Box		
Ultimate 60mm Endoscopic Articulating Reloading Unit	PS9ECAR60N	Black	Very Thick	55mm	6	90	5.0mm	2.35mm	6 / Box		

**Purple Surgical** is a leading independent British manufacturer of high quality, cost-effective single-use surgical instrumentation and devices, partnering with global healthcare providers for positive patient outcomes.

# purple surgical

Ultimate Endoscopic Linear Cutter Staplers with Universal Handles



FROM THE UK'S LEADING MANUFACTURER OF LAPAROSCOPIC INSTRUMENTS AND DEVICES

#### **Contact Us**

Tel: +44 (0)1923 839333 Email: mail@purplesurgical.com www.purplesurgical.com



Great British Manufacturing

Purple Surgical UK Ltd

Purple Surgical International Ltd

2 Chestnut House, Farm Close, Shenley, Herts. WD7 9AD, England

## The Purple Surgical **Ultimate Endoscopic Linear Cutter Staplers**

The Purple Surgical *Ultimate* Endoscopic Linear Cutter Stapler places two, triple-staggered rows of titanium staples, and simultaneously cuts and divides the tissue between the two rows.

**Articulating Reloading Units** 

- The Universal Handle is available in 3 shaft lengths of 6cm, 16cm and 26cm, accommodating 30mm, 45mm and 60mm articulating reloading units.
- The reloading units are available in the following sizes, for varying tissue thicknesses: White (2.5mm)

Blue (3.5mm) Gold (4.2mm) Green (4.8mm) Black (5.0mm)


- · White, Blue and Gold reloading units are compatible with Purple Surgical's **Ultimate** 12mm Trocars.
  - · Green and Black reloading units fit through 15mm Trocars.
  - Each reloading unit houses a sharp knife blade to ensure a clean cut on every firing.
  - · Staple gap control allows for consistent staple formation as the cutting knife proceeds through the transecting tissue.

### Features & Benefits:

- Gentle tissue compression before firing supports repositioning and minimises tissue trauma.
- The cutting knife housed within the reloading unit guarantees an even distance between the anvil and cartridge as it passes through the tissue for consistent, leak free staple formation.
- The cartridge lock-out safety mechanism prevents the reuse of fired reloading units.
- The Universal Handle is available in three size options and is compatible with all Purple Surgical reloading units.
- The Universal Handle can be fired up to 25 times during a single procedure.

• Suitable for: Bariatic Surgery Colorectal Surgery Gynaecological Surgery Paediatric Surgery Thoracic Surgery



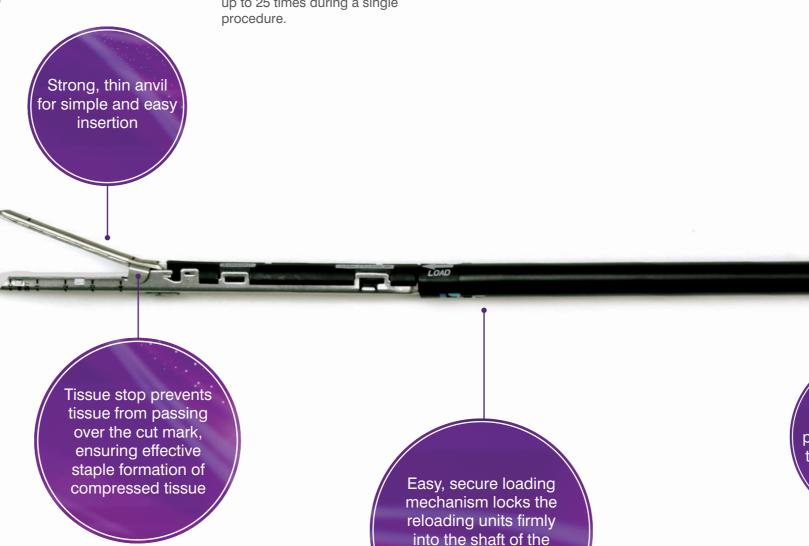

Universal Handle by twisting 45° clockwise

Safety interlock button prevents premature firing and the re-firing of used reloading units

Ergonomic, firing handle with non-slip grip for surgeon comfort and control

Sequential squeezing of the handle is required to fully fire a reloading unit

Single-handed

finger rotation

wheel provides

360° revolutions

Size of No. of Squeezes Reload Unit Required 30mm 45mm 60mm

releases anvil

to allow for the

re-positioning

of tissue

45 30

Gentle closing of anvil to cartridge reduces tissue trauma

> 6 rows of staggered staples